




Wöhlerstr. 2-6
D-55120 Mainz

Tel. 06131/68 66 19
Fax 06131/68 66 20

USt. Id.-Nr.
DE 149057837

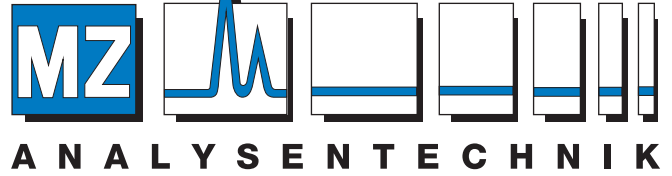

HPLC-Column Quality Certificate

Column: 250 x 4.6 mm [Length x ID]

Part-No.: MZ0831-250046

Material: PerfectSil Target ODS-3 HD

Ser.-No.: 23142939

Particle size: 5 μm

Batch-No.: RF0110

Chromatographic Conditions

Eluent: acetonitrile / water 80/20


Flow rate: 1.00 ml / min

Back pressure: 6.5 MPa

Sample Volume: 0.5 μl

Detection: UV / 254 nm

Temperature: ambient

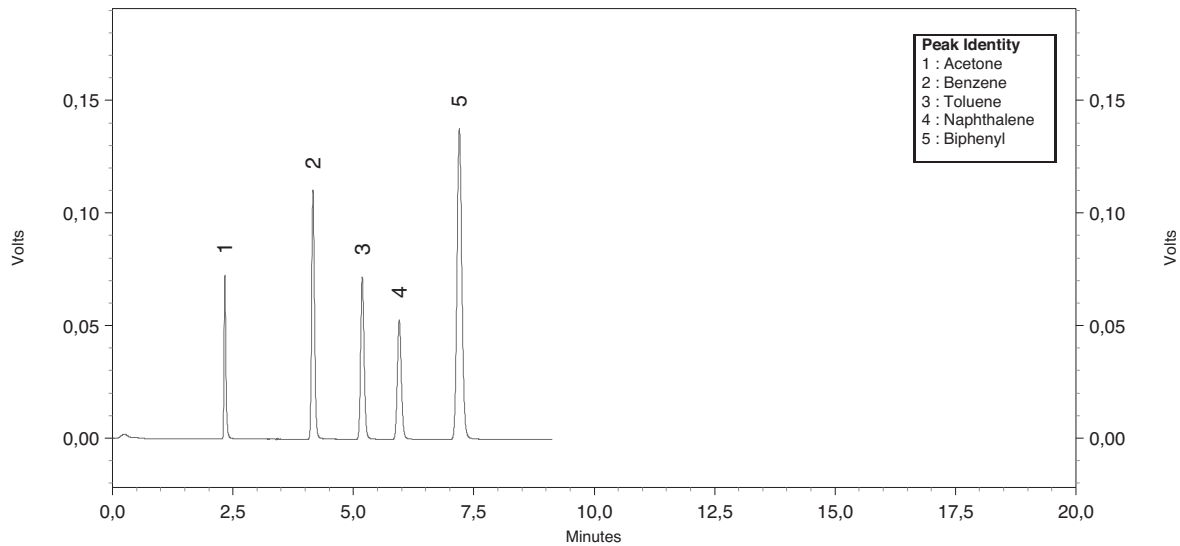

Results

Peak Name	Retention Time [min]	Capacity Factor k'	Theor. Plates / Column	Theor. Plates / m	Asymmetry
Biphenyl	7,20	2,09	26561	106244	1,06

Operator: M. Minz

